



V For *in vitro* Veterinary
Diagnostics only.

Kylt[®]

Kylt[®] IBDV Screening

Real-Time RT-PCR Detection

www.kylt.eu



Kylt® IBDV Screening

Real-Time RT-PCR Detection

Revision No.	Amendments
003	Changed Company name and adjusted layout of the test evaluation

A. General

- Kylt® IBDV Screening kits are intended for the specific detection of viral RNA of Infectious Bursal Disease Virus (IBDV). The kits are suitable for the analysis of samples from birds, such as tissues and organs (e.g. Bursa of Fabricii, spleen, lymphoid tissue), swab samples of the aforementioned tissues and organs and samples from cultural processes with the aforementioned sample material.
- The qualitative testing with Kylt® IBDV Screening kits is based on a duplex Real-Time RT-PCR: In one reaction setting, the RNA target sequences for IBDV as well as for the endogenous control (beta-Actin) are reverse transcribed (Reverse Transcription (RT)) and amplified in parallel with respective primer pairs in the Polymerase Chain Reaction (PCR). Amplified target gene fragments are detected via fluorescently labeled probes during the PCR reaction in real-time (Real-Time PCR). The probes specific for detection of amplified IBDV and the endogenous control target genes are labeled with fluorescent dyes FAM and HEX, respectively, and their emitted fluorescence is separately optically measured by the Real-Time PCR thermal cycler. By means of both individual analyses in one reaction vessel per sample and the Negative Control and Positive Control per run the IBDV-specific status of a sample can be evaluated in the end. This way, results can be achieved within a few hours after sample receipt.
- These kits were developed for use by trained laboratory personnel following standardized procedures. This Direction For Use must be followed strictly.

B. Reagents and Materials

- The following Kylt® IBDV Screening kits are available and comprise the following reagents:

Reagent	Colour of Lid	100 Reactions Article No 31044	25 Reactions Article No 31045	Store at
2x RT-qPCR-Mix	○ transparent	4 x 280 µl	1 x 280 µl	≤ -18 °C
Detection-Mix	● violet	4 x lyophilizate (final 150 µl each)	1 x lyophilizate (final 150 µl each)	≤ -18 °C
Positive Control	● red	4 x lyophilizate (final 50 µl each)	2 x lyophilizate (final 50 µl each)	≤ -18 °C
Negative Control	● blue	1 x 1 ml	1 x 1 ml	≤ -18 °C

- After receipt, the components are immediately stored at ≤ -18 °C. Avoid repeated freezing and thawing of all the reagents and keep them thawed as short as possible.
- If occasional processing of few samples only is expected you may prepare appropriate aliquots of reagents before storage at ≤ -18 °C. Prepare aliquots in such a way that freeze-thaw-cycles are reduced to a maximum of three. The Negative Control can alternatively be stored at +2°C to +8°C.
- The components are to be used within the indicated shelf life (see box label). The components of different batches may not be mixed.
- Before its first use, rehydrate the Positive Control: add 50 µl of Negative Control per vial, briefly incubate at room temperature and mix thoroughly by repeated vortexing. It is recommended to generate aliquots of suitable volumes and store them at ≤ -18 °C.
- The Detection-Mix needs to be stored protected from abundant light. Do not expose to direct (sun)light. Before first use, rehydrate the lyophilized Detection-Mix: add 150 µl of the Negative Control per vial of Detection-Mix, briefly incubate at room temperature and mix by pulse-vortexing. Generate aliquots of suitable volumes and store them at ≤ -18 °C.

C. Equipment and Reagents not included

- This detection method can be used on all commercially available Real-Time PCR thermal cyclers that detect the emitted fluorescence of the fluorescent dyes FAM and HEX (emission 520 and 550 nm, respectively). Note that default normalization option against ROX (e.g. using ABI cyclers) must be deactivated.
- Apart from the disposables, the following further devices are needed and are not included in the Kylt® IBDV Screening kits:
 - RNA preparation kit / protocol (e.g. Kylt® RNA / DNA Purification products)
 - Table top microcentrifuge
 - Vortex
 - Micropipettes covering volumes of 1 µl to 1000 µl
 - Centrifuge for PCR tubes or plates

- Accessory Kylt® products: see chapter F “Related and Accessory Products”.
- We recommend the exclusive use of certified Nuclease-free disposables as well as powder-free protective gloves. Please wear gloves during the entire experimental procedure. Gloves need to be changed frequently, especially after spillage or suspected contaminations.

D. Control Reactions

- The Positive Control allows for control of the specificity and efficiency of the reagents and the reaction itself, including the performance of RT and Real-Time PCR and of the Real-Time PCR thermal cycler.
- The Negative Control allows for exclusion of contaminations. The sample testing is only valid if both, Positive and Negative Controls, are used and verified for validity in every Real-Time PCR run.
- The Internal Control is based on detection of beta-Actin RNA, which is ubiquitous in the cells of the host that the sample is derived from. The beta-Actin RNA/endogenous RNA target gene is co-amplified (channel HEX) with every single reaction and allows for evaluation of sufficient sampling, sample storage and shipment, sample preparation and the Real-Time RT-PCR run itself.
- It is recommended to run one or more of a RNA Isolation Control (RIC) per set of RNA preparation, depending on the total number of samples processed at once. The RIC is a “mock sample” composed of the plain sterile buffer used for raw sample processing. It is randomly placed between the samples, processed like a normal sample and allows to detect potential contaminations of the reagents used (additionally to the Negative Control reaction) as well as for the detection of potential carryover contaminations between individual samples, e.g. during the RNA preparation process.

E. Protocol *(see also „Protocol At A Glance“ at the end of this Direction For Use)*

- The overall protocol of the analysis consists of the following main workflow:
 1. Sample Preparation
 2. RNA Preparation
 3. Reaction Setup, Reverse Transcription and Amplification (Real-Time RT-PCR)
 4. Data Analysis – Validity and Qualitative Result
- We recommend proceeding through the protocol without interruption to avoid potential degradation of the processed samples and reagents. If necessary, you may store the final RNA preparation at $\leq -18\text{ °C}$ or $\leq -70\text{ °C}$ until further processing. Avoid repeated freezing and thawing of the RNA preparations.

1. Sample Preparation

- We recommend pooling of at most five samples or samples from five individuals, respectively, per RNA preparation.
- Pool swabs in a sufficient volume of sterile buffer (e.g. 1 ml of Normal Saline or 0.1 x TE), let the swabs soak for an adequate period of time and finally wash out the swabs by thorough pulse-vortexing. The washed out supernatant is used for RNA preparation.
- Tissue and organ samples are homogenized thoroughly in sterile buffer (see above) and a suitable volume is used for the RNA preparation.
- Material derived from cultural processes, i.e. cell culture supernatant or allantoic fluid, can be used directly for RNA preparation.

2. RNA Preparation

a) Kylt® RNA Preparation (requires Kylt® RNA / DNA Purification products (available separately))

- For detailed information on the RNA preparation process please refer to the manual of Kylt® RNA / DNA Purification products.

b) RNA Preparation by other Methods

- Other kits or in-house methods to purify RNA may be used, as long as the quality leads to satisfactory amplification and detection of the Internal Control RNA.

3. Reaction Setup and Amplification (Real-Time RT-PCR)

- Before each use, briefly vortex and spin down the 2x RT-qPCR-Mix, rehydrated Detection-Mix and Negative Control.
- To determine the total number of reactions needed, count the number of samples and add two more for the Negative Control and the Positive Control (and RIC(s), if processed).
- Prepare the Master-Mix using the components listed below. A larger volume of a ready to use Master-Mix can be prepared and stored at ≤ -18 °C for convenient use over a longer period of time up to the expiry date given on the label. In case of frozen storage the Master-Mix should be aliquoted in such a way that freeze-thaw-cycles are reduced to a maximum of three.
- Vortex, spin down and add 16 μ l of the finalized Master-Mix to each of the PCR tubes or plate wells ("cavities").

Reagent	Volume (μ l)	
	per Reaction	e.g. n=7
2x RT-qPCR-Mix	10 μ l	70 μ l
Detection-Mix	6 μ l	42 μ l
Total Master-Mix	16 μl	112 μl dispense 16 μ l per reaction
RNA (Negative Control / sample RNA / RIC(s) / Positive Control)	4.0 μ l	
Total Reaction	20.0 μl	

- Keep exposure of the 2x RT-qPCR-Mix, Detection-Mix and prepared Master-Mix to (sun)light as short as possible and return it back to appropriate storage temperature right after application. Avoid the formation of bubbles when pipetting the Master-Mix, samples and controls.
- Add 4 µl of the Negative Control to the corresponding cavity and seal it individually, if possible.
- Add 4 µl of each RNA preparation (including RIC(s), if processed) to the corresponding cavities and seal them individually, if possible.
- To minimize risk of potential cross-contaminations, 4 µl of the Positive Control are added to the corresponding cavity after all previous samples and control reactions are set up. Before each use, briefly vortex and spin down the rehydrated Positive Control (see also chapter B “Reagents and Materials”).
- If not already done, finally seal the cavities. It is recommended to briefly spin them down before the start of the Real-Time PCR run.
- Place the cavities in the Real-Time PCR thermal cycler and run the test with Kylyt® Profile I or III as given below.

Kylyt® Profile I or III				
Step No	Description	Temperature	Duration	
1	Reverse Transcription	50 °C	10 min	
2	Activation of Polymerase	95 °C	1 min	
3	Denaturation	95 °C	10 sec	} 42 cycles
4	Annealing & Extension	60 or 55 °C	1 min	
5	Fluorescence Detection	channels FAM and HEX		

- Kylyt® Profile I allows for combined run of this and most other Kylyt® RT-qPCR detection methods as well as Kylyt® PCR detection products.
- Kylyt® Profile III allows for combined run of the Kylyt® RT-qPCR detection products for IBV variants and IBDV.
- In the event of a combined Real-Time RT-PCR run, make sure all necessary channels are detected.
- Please follow the specified instructions of your Real-Time PCR thermal cycler as recommended by the manufacturer.

4. Data Analysis – Validity and Qualitative Result

General

- The amplification data can be processed automatically using the specific software tool of your Real-Time PCR thermal cycler. Alternatively, the threshold can be set manually considering the following directions: The threshold should cross the FAM-curve and the HEX-curve in the linear increase of their slope (log scaling of the y-axis). By setting the threshold, the crossing points with the HEX- and FAM-curves determine the respective cycle threshold (Ct), which is negatively correlated with the initial concentration of copies of the target genes in the Real-Time RT-PCR reaction.
- Only curves with the typical exponential amplification, meaning the curve of the raw data shows a flat baseline at the beginning, followed by a clear (exponential) slope in fluorescence and possibly reaching a plateau-phase (y-axis in log scaling), should be regarded as positive.

- The actual test analysis starts with the validity check of the entire Real-Time RT-PCR run. Afterwards, by means of the Internal Control the validity of each sample reaction and its true test result can be verified according to the Ct-value of the Internal Control channel (HEX). Finally, the IBDV-specific status of each sample is analyzed (FAM).

Test Evaluation - Control Reactions

- The **Real-Time RT-PCR test run** is only **valid** if the curves of the control reactions can be evaluated as follows:

Control Reactions	Channel			
	HEX (beta-actin)		FAM (IBDV)	
Negative Control	negative	Ct > 35	negative	Ct > 42
Positive Control	positive	Ct ≤ 35	positive	Ct > 15 and ≤ 35

- If one or more of a RNA Isolation Control (RIC(s)) is processed, its FAM- and HEX-curves must be negative (Ct > 35 for HEX-curves).

Test Evaluation - Samples

Target	Channel	Signal		
Internal Control (beta-actin)	HEX	positive	positive / negative	negative
IBDV	FAM	negative	positive	negative
The sample is IBDV		negative	positive	inhibited

- A **sample** is **negative for Infectious Bursal Disease Virus**, if its HEX-curve is positive (Ct ≤ 40), but its FAM-curve is negative.
- A **sample** is **positive for Infectious Bursal Disease Virus**, if its FAM-curve is positive (Ct ≤ 42), independent of the HEX-curve.
- A **sample** is **inhibited** if neither the FAM-curve nor the HEX-curve are positive.
- **Recommendation:** In the case of an inhibited sample you may repeat the test by using e.g. 1:4 dilution of the respective RNA preparation. The Negative Control is used as the diluting agent. Preferably, the entire RNA preparation process is repeated using Kylt® RNA/DNA Purification products or appropriate alternative.
- Convenient and reliable sample data entry, Real-Time PCR start, final qualitative analysis and documentation can be conducted with the Kylt® Software, please inquire.

F. Related and Accessory Products

Product	Article No	Content	Description
Kylt® IBDV genogroup A1 (classical)	31667/ 31668	100 / 25	PCR Kit for the specific detection of viral RNA of <i>Infectious Bursal Disease Virus (IBDV)</i> Genogroup A1.
Kylt® IBDV genogroup A3 (very virulent)	31601/ 31602	100 / 25	PCR Kit for the specific detection of viral RNA of <i>Infectious Bursal Disease Virus (IBDV)</i> Genogroup A3.
Kylt® IBDV genogroup A7 (early Australian)	31669/ 31670	100 / 25	PCR Kit for the specific detection of viral RNA of <i>Infectious Bursal Disease Virus (IBDV)</i> Genogroup A7.
Kylt® RNA / DNA Purification	31315	50	Combined RNA and DNA purification from veterinary samples (spin-column based).
Kylt® RNA / DNA Purification HTP	31826	4 x 96	Magnetic bead based combined RNA and DNA purification kit for veterinary diagnostic samples. Suitable for Kylt® Purifier and Kylt® Purifier 48.
Kylt® Purifier	31436	1 unit	Purification system for magnetic bead based kits. Up to 96 samples are processed in under 30 minutes. Intended for high-throughput laboratories.
Kylt® Purifier 48	31748	1 unit	Purification system for magnetic bead based kits. Up to 48 samples are processed in under 30 minutes. Intended for low to medium throughput laboratories.
Kylt® Purifier Spin Tips	31434	5 Sets	Plate with 96 separate spin tips, used by the Kylt® Purifier to mix the well contents by stirring. Sufficient for 480 samples.
Kylt® Purifier Plates	31435	20 Plates	Plates to be used for the several reactions and reagents during automated nucleic acid purification. Sufficient for 320 to 480 samples (depending on device and protocol).

G. Ordering information

For a fast and efficient service please send your order to orders.kylt-de@san-group.com and please provide the following information:

- Delivery address
- Invoice address
- Purchaser contact telephone number
- End user name and telephone number (if different)
- Purchase order number
- Product name and catalogue number
- Quantity and size of products
- Indicate if your account is VAT exempt

Production:

SAN Group Biotech Germany GmbH | Muehlenstr. 13 | 49685 Hoeltinghausen | Germany
www.kylt.eu | kylt-de@san-group.com

Development, manufacturing and distribution of Kylt® *In-Vitro* Diagnostica is certified according to ISO 9001:2015.

Kylt® is a registered trademark.

For veterinary use only. For *in vitro* use only. Regulatory requirements vary by country, not all of the products described herein may be available in your geographic area.

© 2023 SAN Group Biotech Germany GmbH. All rights reserved. The trademark mentioned herein is the property of SAN Group Biotech Germany GmbH or their respective owners.




PROTOCOL AT A GLANCE

Real-Time RT-PCR Setup

1


Prepare a Master-Mix*

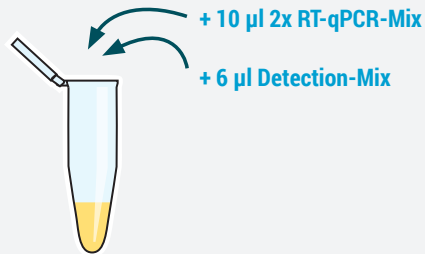

* please refer to chapter 3


2

Pulse-vortex and spin down

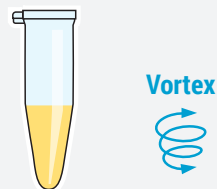


3

Dispense Master-Mix and add 4 µl NC, RNA preparation, PC

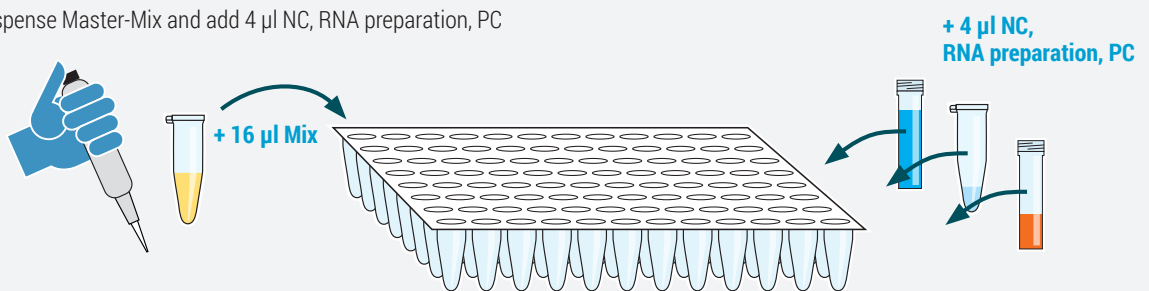

4

Seal cavities, spin down (recommended), and start cycler

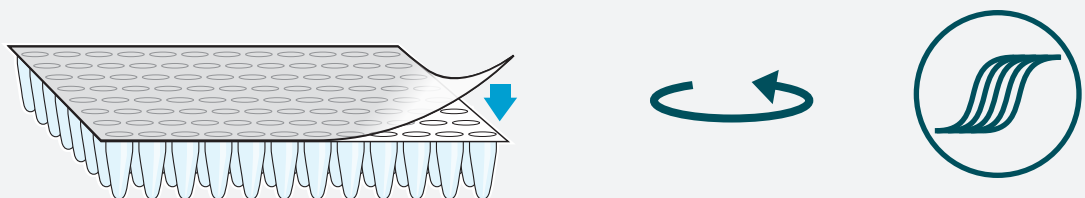

5

Analysis

